Elabscience®

Annexin V-EGFP Azid-Free Lyophilized Powder

Lyophilized Powder

Cat. No: E-CK-A119U		Size: 25 µg	Size: 25 µg/50 µg/100 µg/200 µg				
Cat.	Products	25 μg	50 µg	100 µg	200 µg	Storage	
E-CK-A119U	Annexin V-EGFP Azide-Free	25 μg	25 μg×2	100 µg	100 μg×2	-20 °C, shading light	


Storage

The lyophilized powder is stable for at least two years when stored at -20 $^{\circ}$ C in the dark after arrival.

Introduction

Elabscience[®] Annexin V-EGFP Azide-Free Lyophilized Powder is developed for long-term storage and userdefined concentrations of EGFP-conjugated Annexin V. It contains no preservatives like azides or Proclin. Annexin V is a member of the annexin family, which binds to phosphatidylserine (PS) in a calcium-dependent manner. The EGFP-conjugated format of this protein, Annexin V-EGFP, binds specifically to the PS displayed on the outer leaflet of apoptotic cell membrane. The signal of Annexin V-EGFP can be collected with flow cytometry or fluorescence microscopy to detect cell apoptosis.

Jurkat cells were treated with 5 μ M Camptothecin and detected with this reagent and PI.

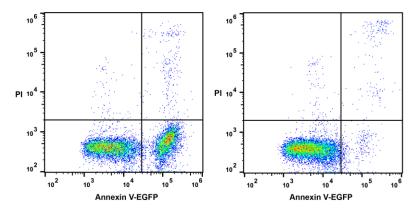

Jurkat cells were cultured with (Left) or without (Right) 5 µM Camptothecin for 4 h. Annexin V-EGFP single-positive cells were early apoptotic cells, Annexin V-EGFP and PI doublepositive cells were necrotic or late apoptotic cells, and PI single-positive cells were naked nuclei.

Reconstitution

Reconstitute the reagent with sterile deionized water. Make sure the final concentration is less than 1 mg/mL. (This concentration is high enough for stock solution. In our internal tests, working solutions with $1\sim2$ µg/mL Annexin V-EGFP have best results for flow cytometry.) For example, at least 10 µL water should be added to 10 µg Annexin V-EGFP. Otherwise, the lyophilized cake may not be thoroughly soaked in water.

Cautions

- 1. This kit is for reseach use only.
- 2. Detect apoptosis as soon as possible after staining to avoid the increase number of apoptosis or necrosis. Avoid extended exposure of the samples to direct light to protect the fluorophores from quenching.
- 3. For your safety and health, please wear the lab coat and disposable gloves before the experiments.

For Research Use Only

Toll-free: 1-888-852-8623 Web: <u>www.elabscience.com</u>