

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

# Elab Fluor® 488 Anti-Mouse CD71 Antibody[R17 217.1.3/TIB-219]

Catalog No.E-AB-F1093LReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

## **Antigen Information**

Alternate Names Transferrin receptor protein 1,Tfrc,TR,TfR,TfR1,Trfr,CD71

**Uniprot ID** Q6235

**Background** CD71 is a 95 kD type II heterodimeric transmembrane glycoprotein that is also known as T9 and

transferrin receptor. CD71 is expressed on proliferating cells, reticulocytes, and erythroid precursors. Its expression is very low on resting leukocytes. CD71 plays a role in the control of cellular proliferation by facilitating the uptake of iron via ferrotransferrin binding and the

recycling of apotransferrin to the cell surface.

### **Product Details**

Form Liquid

Size 50Tests/100Tests/100Tests×2

**Clone No.** R17 217.1.3/TIB-219

**Host** Rat

IsotypeRat IgG2a, κReactivityMouseApplicationFCM

Isotype Control Elab Fluor<sup>®</sup> 488 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09832L]

**Storage Buffer** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

## **Fluorophore**


Conjugation: Elab Fluor® 488

Elab Fluor® 488 is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 520 nm (e.g., a 525/40 nm bandpass filter).

## Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μL of antibody per test (million cells in 100 μL staining volume or per 100 μL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

#### **Product data**

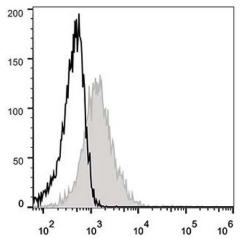

C57BL/6 murine splenocytes are stained with Elab Fluor<sup>®</sup> 488 Anti-Mouse CD71 Antibody (filled gray histogram). Unstained splenocytes (empty black histogram) are used as control.

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/Listdetail-459742.html