

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

# Elab Fluor® 647 Anti-Mouse MHC II (I-A/I-E) Antibody[M5/114]

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

# **Antigen Information**

Alternate Names H2-Ab1/Eb1, Major histocompatibility protein class II beta chain, MHC class II H2-IA-beta-psi, I-

E beta MHC class II, MHC class II

**Uniprot ID** P14483,O78196

**Background** These class II molecules are expressed on antigen presenting cells (including B cells) and a subset

of T cells from H-2b,d,q,r bearing mice and are involved in antigen presentation to T cells

expressing CD3/TCR and CD4 proteins.

#### **Product Details**

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. M5/114 Host Rat

IsotypeRat IgG2b, κReactivityMouseApplicationFCM

**Isotype Control** Elab Fluor<sup>®</sup> 647 Rat IgG2b, κ Isotype Control[LTF-2] [Product E-AB-F09842M]

**Storage Buffer** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

# **Fluorophore**


Conjugation: Elab Fluor® 647

Elab Fluor® 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).

### Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μL of antibody per test (million cells in 100 μL staining volume or per 100 μL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

#### **Product data**

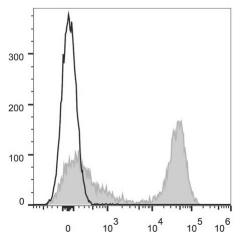

C57BL/6 murine splenocytes are stained with Elab Fluor<sup>®</sup> 647 Anti-Mouse MHC II (I-A/I-E) Antibody (filled gray histogram). Unstained splenocytes (empty black histogram) are used as control.

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5594.html">https://www.elabscience.com/List-detail-5594.html</a>
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/Listdetail-459742.html

Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation