

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

# PerCP Anti-Mouse CD48 Antibody[HM48-1]

Catalog No.E-AB-F1017FReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

## **Antigen Information**

Alternate Names CD48 antigen, Cd48, BCM1 surface antigen, BLAST-1, HM48-1, MRC OX-45 surface

antigen, SLAMF2, sgp-60, CD48

Uniprot ID P18181

**Background** CD48 is a 45 kD GPI-anchored glycoprotein also known as BCM1, Blast-1 (human), and OX-45

(rat). It is a member of the Ig superfamily, expressed on T and B cells and

monocytes/macrophages. It plays a role in adhesion and T cell recognition. The primary ligands

for CD48 are CD2 and CD244.

### **Product Details**

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. HM48-1

Host Armenian Hamster
Isotype Armenian Hamster IgG

**Reactivity** Mouse **Application** FCM

**Isotype Control** [Product E-AB-F09852F]

**Storage Buffer** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

# **Fluorophore**


Conjugation: PerCP

PerCP is designed to be excited by the blue laser (488 nm) and detected using an optical filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).

## Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 µL staining volume or per 100 µL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

#### **Product data**

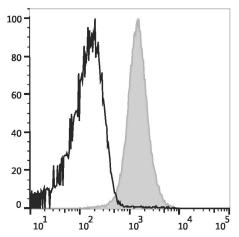

C57BL/6 murine splenocytes are stained with PerCP Anti-Mouse CD48 Antibody (filled gray histogram) or PerCP Armenian Hamster IgG Isotype Control (empty black histogram).

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/Listdetail-459742.html

Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation