# **Elabscience**®

# FITC Anti-Human CD134 Antibody[Ber-ACT35]

Catalog No.E-AB-F1153CStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityHApplicationsF

Human FCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

#### **Antigen Information**

| Alternate Names | TXGP1L,TNFRSF4,ACT35 antigen,OX40L receptor,CD134  |
|-----------------|--|
| Uniprot ID      | P43489   |
| Background      | CD134, also known as OX40 and TNFRSF4, is a 50 kD type I transmembrane glycoprotein. It is a member of the TNF receptor family. OX40 is expressed on activated T lymphocytes including Th1, Th2, Th17, and Treg cells. The interaction of OX40 with OX40L results in B cell proliferation and antibody secretion, regulation of primary T cell expansion, and T cell survival. OX40 influences the size of the T cell memory pool and regulation of CD4+ T cell tolerance. |

## **Product Details**

| Form                | Liquid  |
|---------------------|---|
| Size                | 20Tests/100Tests/100Tests×2   |
| Clone No.           | Ber-ACT35   |
| Host                | Mouse   |
| Isotype             | Mouse IgG1, κ   |
| Reactivity          | Human   |
| Application         | FCM   |
| Isotype Control     | FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C]                          |
| Storage Buffer      | Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. |
| Shipping            | Biological ice pack at 4 °C   |
| Stability & Storage | Keep as concentrated solution.  |
|                     | Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.                |
|                     | This product is guaranteed up to one year from purchase.                                    |

For Research Use Only

Thank you for your recent purchase. If you would like to learn more about antibodies, please visit www.elabscience.com. Focus on your research Service for life science

Applications:Activ- Activation; Block- Blocking; Separation- Cell Separation ; Cell Sep-Neg- Cell Separation by Negative Selection; FA-Functional Assay; Neut- Neutralization; Stim- Stimulation; FCM- Flow Cytometry; ICFCM: Intracellular Staining for Flow Cytometry; WB-Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation

# **Elabscience**®

# Fluorophore


#### **Conjugation:** FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

## **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

# **Product data**

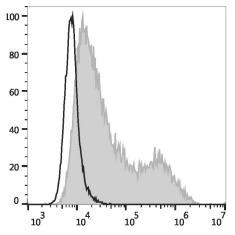

CHO cells transiently transfected with pcDNA3.1 plasmid encoding Human CD134 gene are stained with FITC Anti-Human CD134 Antibody (filled gray histogram) or FITC Mouse IgG1,  $\kappa$  lsotype Control (empty black histogram).

## **Related Information**

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>

For Research Use Only

Thank you for your recent purchase. If you would like to learn more about antibodies, please visit www.elabscience.com. Focus on your research Service for life science

Applications:Activ- Activation; Block- Blocking; Separation- Cell Separation ; Cell Sep-Neg- Cell Separation by Negative Selection; FA-Functional Assay; Neut- Neutralization; Stim- Stimulation; FCM- Flow Cytometry; ICFCM: Intracellular Staining for Flow Cytometry; WB-Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation