Elabscience®

FITC Anti-Human IFN-γ Antibody[B27]

Catalog No.E-AB-F1196CStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityHumanApplicationsICFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names	Interferon gamma,IFNG,IFN-gamma,IFNγ
Uniprot ID	P01579
Gene ID	3458
Background	Interferon- γ is a potent multifunctional cytokine which is secreted primarily by activated NK cells and T cells. Originally characterized based on anti-viral activities, IFN- γ also exerts anti- proliferative, immunoregulatory, and proinflammatory activities. IFN- γ can upregulate MHC class I and II antigen expression by antigen-presenting cells. The B27 antibody reacts with the human interferon- γ . The B27 antibody can neutralize the bioactivity of natural or recombinant IFN- γ .

Product Details

Form	Liquid
-	1
Size	20Tests/100Tests/100Tests×2
Clone No.	B27
Host	Mouse
Isotype	Mouse IgG1, ĸ
Reactivity	Human
Application	ICFCM
Isotype Control	FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

For Research Use Only

Thank you for your recent purchase. If you would like to learn more about antibodies, please visit www.elabscience.com. Focus on your research Service for life science

Applications:Activ- Activation; Block- Blocking; Separation- Cell Separation ; Cell Sep-Neg- Cell Separation by Negative Selection; FA-Functional Assay; Neut- Neutralization; Stim- Stimulation; FCM- Flow Cytometry; ICFCM: Intracellular Staining for Flow Cytometry; WB-Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation

Elabscience®

Fluorophore


Conjugation: FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

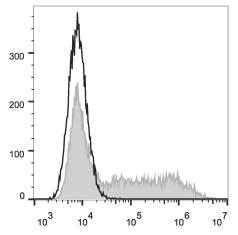

HEK293T cells transiently transfected with pcDNA3.1 plasmid encoding Human IFN- γ gene are stained with FITC Anti-Human IFN- γ Antibody (filled gray histogram) or FITC Mouse IgG1, κ lsotype Control (empty black histogram).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Intracellular Antigens for Flow Cytometry https://www.elabscience.com/List-detail-5570.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>

For Research Use Only Thank you for your recent purchase

If you would like to learn more about antibodies, please visit www.elabscience.com.

Focus on your research Service for life science

Applications:Activ- Activation; Block- Blocking; Separation- Cell Separation ; Cell Sep-Neg- Cell Separation by Negative Selection; FA-Functional Assay; Neut- Neutralization; Stim- Stimulation; FCM- Flow Cytometry; ICFCM: Intracellular Staining for Flow Cytometry; WB-Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation