

Elab Fluor® 647 Anti-Human CD74 Antibody[LN2]

Catalog No.	E-AB-F1072M	Reactivity	Human
Storage	Store at 2~8°C, Avoid freeze / thaw cycles	Applications	FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names	HLA class II histocompatibility antigen gamma chain,CD74,Ia antigen-associated invariant chain,Ii,p33,DHLAG
Uniprot ID	P04233
Background	CD74 is a type II transmembrane glycoprotein also known as MHC class II associated invariant chain, invariant chain, Ii, MHC class II chaperone, and MIF receptor. CD74 exists in four isoforms with molecular masses of 33, 35, 41, and 43 kD, depending on genetic splicing. CD74 is primarily expressed on antigen presenting cells, including B cells, monocytes/macrophages, dendritic cells, and Langerhans cells. It is also expressed by activated T cells and activated endothelial and epithelial cells as well as carcinomas of lung, renal, gastric and thymic origin. The primary function of CD74 is intracellular sorting of MHC class II molecules and regulation of exogenous peptide loading onto MHC class II. It is also involved in the modulation of B cell differentiation and positive selection of CD4+ T cells. It has been reported that CD74 binds MIF (macrophage migration inhibitory factor) and signals through CD44 to regulate innate and adaptive immunity. It is also reported that H. pylori infection occurs through urease B binding of CD74 on gastric epithelial cells, inducing gastric epithelial cell apoptosis, NF-κB activation, and IL-8 production.

Product Details

Form	Liquid
Size	20Tests/100Tests/100Tests×2
Clone No.	LN2
Host	Mouse
Isotype	Mouse IgG1, κ
Reactivity	Human
Application	FCM
Isotype Control	Elab Fluor® 647 Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792M]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light.Do not freeze. This product is guaranteed up to one year from purchase.

For Research Use Only

Thank you for your recent purchase.

If you would like to learn more about antibodies,please visit www.elabscience.com.

**Focus on your research
Service for life science**

Applications:Activ- Activation; Block- Blocking; Separation- Cell Separation ; Cell Sep-Neg- Cell Separation by Negative Selection; FA- Functional Assay; Neut- Neutralization; Stim- Stimulation; FCM- Flow Cytometry; ICFM: Intracellular Staining for Flow Cytometry; WB- Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation

Fluorophore


Conjugation: Elab Fluor® 647

Elab Fluor® 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 µL staining volume or per 100 µL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

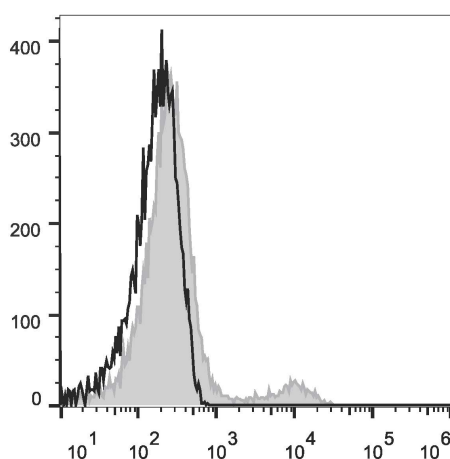

Human peripheral blood lymphocytes are stained with Elab Fluor® 647 Anti-Human CD74 Antibody (filled gray histogram). Unstained lymphocytes (empty black histogram) are used as control.

Related Information

1. Sample Preparation for Flow Cytometry <https://www.elabscience.com/List-detail-5594.html>
2. Staining Cell Surface Targets for Flow Cytometry <https://www.elabscience.com/List-detail-5568.html>
3. Flow Cytometry Troubleshooting Tips <https://www.elabscience.com/List-detail-5593.html>
4. How to select the appropriate detection channel through the spectrogram? <https://www.elabscience.com/List-detail-459742.html>

For Research Use Only

Thank you for your recent purchase.

If you would like to learn more about antibodies, please visit www.elabscience.com.

**Focus on your research
Service for life science**

Applications: Activ- Activation; Block- Blocking; Separation- Cell Separation ; Cell Sep-Neg- Cell Separation by Negative Selection; FA- Functional Assay; Neut- Neutralization; Stim- Stimulation; FCM- Flow Cytometry; ICFM: Intracellular Staining for Flow Cytometry; WB- Western Blotting; IHC- Immunohistochemistry; IF- Immunofluorescence; IP- Immunoprecipitation