

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

Elab Fluor® 647 Anti-Mouse IgD Antibody[11-26c.2a]

Catalog No.E-AB-F1189UMReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Immunoglobulin heavy constant delta,IGHD,Ig delta chain C region

Uniprot ID P01881

Background Surface IgD is an important B cell differentiation marker.

Product Details

 Form
 Liquid

 Concentration
 0.5 mg/mL

 Size
 25μg/100μg

 Clone No.
 11-26c.2a

 Host
 Rat

IsotypeRat IgG2a, κReactivityMouseApplicationFCM

Isotype Control Elab Fluor[®] 647 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09833M]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

Fluorophore


Conjugation: Elab Fluor® 647

Elab Fluor $^{\odot}$ 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the reagent to obtain optimal results [The recommended concentration is 0.1-1 $\mu g/10^6$ cells in $100~\mu L$ volume].

Product data

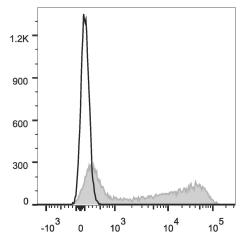

C57BL/6 murine splenocytes are stained with Elab Fluor[®] 647 Anti-Mouse IgD Antibody (filled gray histogram).

Unstained splenocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html