

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

PerCP/Cyanine5.5 Anti-Mouse CD279/PD-1 Antibody[29F.1A12]

Catalog No. E-AB-F1131J **Reactivity** Mouse **Storage** Store at 2~8°C, Avoid freeze / thaw cycles **Applications** FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names PD-1, Programmed Death-1

Uniprot ID Q02242

Background CD279, also known as programmed death-1 (PD-1), is a 50-55 kD glycoprotein belonging to the

CD28 family of the Ig superfamily. PD-1 is expressed on activated splenic T and B cells and thymocytes. It is induced on activated myeloid cells as well. PD-1 is involved in lymphocyte clonal selection and peripheral tolerance through binding its ligands, B7-H1 (PD-L1) and B7-DC (PD-L2). It has been reported that PD-1 and PD-L1 interactions are critical to positive selection and play a role in shaping the T cell repertoire. PD-L1 negative costimulation is essential for

prolonged survival of intratesticular islet allografts.

Product Details

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. 29F.1A12 Host Rat

IsotypeRat IgG2a, κReactivityMouseApplicationFCM

Isotype Control PerCP/Cvanine5.5 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09832J]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

Fluorophore


Conjugation: PerCP/Cyanine5.5

PerCP/Cyanine 5.5 is designed to be excited by the blue laser (488 nm) and detected using an optical filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μL of antibody per test (million cells in 100 μL staining volume or per 100 μL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

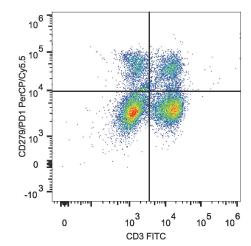

C57BL/6 murine splenocytes are stained with PerCP/Cyanine5.5 Anti-Mouse CD279/PD-1 Antibody FITC Anti-Mouse CD3 Antibody.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/Listdetail-459742.html