

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

PerCP/Cyanine5.5 Anti-Human CD9 Antibody[HI9a]

Catalog No.E-AB-F1086JReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names CD9 antigen, CD9,5H9 antigen, Cell growth-inhibiting gene 2 protein, Leukocyte antigen

MIC3,MRP-1,Tspan-29,p24,CD9,MIC3, TSPAN29

Uniprot ID P21926

Background CD9 is a 24 kD type III transmembrane protein also known as tetraspanin, MRP-1 and DRAP-24.

It is a member of the tetraspan family (spanning the membrane four times) found on platelets, B cell progenitors, activated lymphocytes, granulocytes, endothelial cells and epithelial cells. CD9 induces adhesion, platelet aggregation, and B cell development. CD9 has been shown to associate

with CD63, CD81, CD82, and CD36 and to bind to β1 integrins.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. HI9a Host Mouse

Isotype Mouse IgG1, κ

Reactivity Human **Application** FCM

Isotype ControlPerCP/Cyanine5.5 Mouse IgG1. κ Isotype Control[MOPC-21] [Product E-AB-F09792J]Storage BufferPhosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

Fluorophore


Conjugation: PerCP/Cyanine5.5

PerCP/Cyanine5.5 is designed to be excited by the blue laser (488 nm) and detected using an optical filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

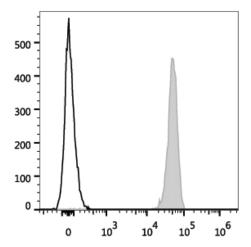

Human platelets are stained with PerCP/Cyanine5.5 Anti-Human CD9 Antibody (filled gray histogram). Unstained platelets (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html