

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

PE Anti-Human CD133 Antibody[W6B3C1]

Catalog No.E-AB-F1268DReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Prominin-1, PROM1, AC133

Uniprot ID O43490

Background CD133, also known as Prominin-1 and AC133 antigen, is a 120 kD pentaspan glycoprotein with 5

transmembrane domains. CD133 was initially described as a surface antigen specific for human hematopoietic stem cells and as a marker for murine neuroepithelial cs and some embryonic epithelia. Later on, CD133 was found on other stem cells, including endothelial progenitor cells, glioblastomas, neuronal, and glial stem cells. In addition to stem cells for normal tissue, CD133 was found on cancer cells, such as some leukemia cells and brain tumor cells. Although the biological function of CD133 is not completely understood, CD133 has been extensively used as a

stem cell marker for normal and cancerous tissues.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No.W6B3C1HostMouseIsotypeMouse IgG1, κReactivityHuman

Application FCM

Isotype Control PE Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792D]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

Fluorophore


Conjugation: PE

PE is designed to be excited by the Blue (488 nm), Green (532 nm) and Yellow-Green (561 nm) lasers and detected using an optical filter centered near 575 nm (e.g., a 585/42 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

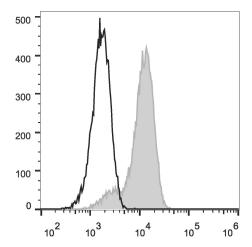

HCT116 cells are stained with PE Anti-Human CD133 Antibody (filled gray histogram). Unstained cells (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Intracellular Antigens for Flow Cytometry https://www.elabscience.com/List-detail-5570.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html