

Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

# APC Anti-Human CD57 Antibody[TB01]

Catalog No.E-AB-F1237EReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

## **Antigen Information**

Alternate Names Galactosylgalactosylxylosylprotein 3-beta-glucuronosyltransferase 1,B3GAT1,GlcAT-P,GlcUAT-

P.GLCATP

**Background** CD57, also known as TB01, NK-1, and Leu-7 is a 100-115 kD oligosaccharide antigenic

determinant expressed on a variety of proteins, lipids, and chondroitin sulfate proteoglycans. CD57 is expressed on a subset of peripheral blood lymphocytes, including NK cells and CD8+ T cells, and is also expressed on neural cells and striated muscle. CD57 is not expressed on red blood cells, granulocytes, monocytes, or platelets. While the function of CD57 is unknown, binding to L-selectin, P-selectin, and a fragment of laminin suggests that CD57 may be involved in cell-matrix interactions. CD57 is increased in some disease states associated with CD4/CD8

imbalances (AIDS, autoimmune disease, viral infections, and allograft transplants).

#### **Product Details**

Form Liquid

Size 20Tests/50Tests/100Tests/200Tests

Clone No. TB01

Host Mouse

Isotype Mouse IgM, κ

Reactivity Human

Application FCM

**Isotype Control**APC Mouse IgM, κ Isotype Control[MM-30] [Product E-AB-F09782E]

Storage Buffer PBS with 0.05% Proclin300, 1% BSA

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com


## **Fluorophore**

Conjugation: APC

### Recommended usage

Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 µL staining volume or per 100 µL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

#### **Product data**

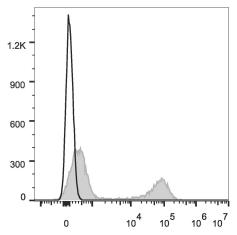

Human peripheral blood lymphocytes are stained with APC Anti-Human CD57 Antibody (filled gray histogram). Unstained lymphocytes (empty black histogram) are used as control.

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5594.html">https://www.elabscience.com/List-detail-5594.html</a>
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/Listdetail-459742.html