


Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

PTGES3 Polyclonal Antibody

Catalog No.E-AB-66063ReactivityHStorageStore at -20°C. Avoid freeze / thaw cycles.HostRabbitApplicationsIHCIsotypeIgG

Note: Centrifuge before opening to ensure complete recovery of vial contents.

Images

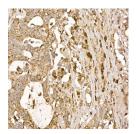

Immunohistochemistry of paraffinembedded human colon carcinoma using P23/PTGES3 Polyclonal Antibody at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.



Immunohistochemistry of paraffinembedded human colon using P23/PTGES3 Polyclonal Antibody at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.



For Research Use Only

Thank you for your recent purchase

If you would like to learn more about antibodies, please visit www.elabscience.com.

Immunogen Information

Immunogen Recombinant fusion protein of human PTGES3

GeneID 10728 Swissprot Q15185

Synonyms PTGES3,P23,TEBP,cPGES

Product Information

Calculated MW 14kDa/16kDa/18kDa

Buffer PBS with 0.02% sodium azide,50% glycerol,pH7.3.

Purify Affinity purification
Dilution IHC 1:50-1:200

Background

This gene encodes an enzyme that converts prostaglandin endoperoxide H2 (PGH2) to prostaglandin E2 (PGE2). This protein functions as a co-chaperone with heat shock protein 90 (HSP90), localizing to response elements in DNA and disrupting transcriptional activation complexes. Alternative splicing results in multiple transcript variants. There are multiple pseudogenes of this gene on several different chromosomes.

Focus on your research Service for life science



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

Immunohistochemistry of paraffinembedded mouse testis using P23/PTGES3 Polyclonal antibody at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Focus on your research Service for life science