


Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

FMR1 Polyclonal Antibody


Catalog No.E-AB-61039ReactivityH,M,RStorageStore at -20°C. Avoid freeze / thaw cycles.HostRabbitApplicationsWB,IHCIsotypeIgG

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Images

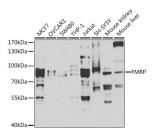

Western blot analysis of extracts of various cell lines using FMRP Polyclonal Antibody at 1:1000 dilution.

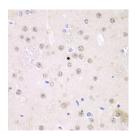

Immunohistochemistry of paraffinembedded rat brain using FMRP Polyclonal Antibody at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

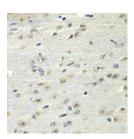

Immunohistochemistry of paraffinembedded mouse brain using FMRP Polyclonal antibody at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS

Immunogen Information

Immunogen Recombinant fusion protein of human FMR1

GeneID 2332 Swissprot Q06787

Synonyms FMRP,FRAXA,POF,POF1,FMR1,POFX

Product Information

Calculated MW 58kDa/61kDa/66kDa/68kDa/69kDa/70kDa/71kDa

Observed MW 85kDa

Buffer PBS with 0.02% sodium azide,50% glycerol,pH7.3.

Purify Affinity purification

Dilution WB 1:500-1:2000,IHC 1:50-1:200

Background

The protein encoded by this gene binds RNA and is associated with polysomes. The encoded protein may be involved in mRNA trafficking from the nucleus to the cytoplasm. A trinucleotide repeat (CGG) in the 5' UTR is normally found at 6-53 copies, but an expansion to 55-230 repeats is the cause of fragile X syndrome. Expansion of the trinucleotide repeat may also cause one form of premature ovarian failure (POF1). Multiple alternatively spliced transcript variants that encode different protein isoforms and which are located in different cellular locations have been described for this gene.

For Research Use Only

Thank you for your recent purchase

If you would like to learn more about antibodies, please visit www.elabscience.com.

Focus on your research Service for life science



Tel:240-252-7368(USA) Fax: 240-252-7376(USA) techsupport@elabscience.com Website: www.elabscience.com

buffer pH 7.2 before commencing with IHC staining protocol.